F21 Maverick™ Xtreme AT

DE Gebrauchsanweisung	3
EN Instructions for use	6
FR Instructions d'utilisation.....	8
NL Gebruiksaanwijzing	11

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Maverick™ Xtreme AT Instructions for Use

Product Number: F21


Assembly

The Maverick Xtreme AT foot module is pre-assembled consisting of glass composite components (keel and sole plate), a Spectra™ sock, and a foot shell. Stiffening bumpers for increasing heel stiffness are provided. After dynamic alignment, torque pyramid adjustment screws to the manufacturer's specifications. Secure pyramid adjustment screws with thread locking adhesive (e.g., Loctite 242).

Bench Alignment

Prior to donning the prosthesis:

- Plantar flex/Dorsiflex the foot to match the shoe heel height (foot designed for shoes with 1 cm heel height).
- Adduct/Abduct the socket to provide appropriate frontal plane angle.
- Flex/Extend the socket to provide appropriate sagittal plane angle.
- Move the socket linearly to ensure the weight line falls along the anterior edge of the pylon (see illustration).

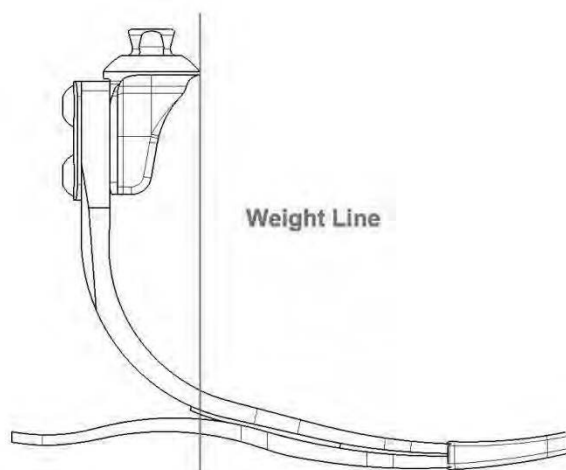Dynamic Alignment

During *loading response*, the heel lever stores energy and releases it during *midstance*. This action provides momentum for the keel to store energy and release it during *terminal stance*. To optimize the heel to toe rollover motion, adjust the following variables:

- Anterior/posterior foot placement
- Dorsiflexion/plantar flexion
- Heel stiffness

Troubleshooting

Heel too soft

Symptoms

- Foot flat occurs too rapidly
- Toe feels excessively stiff
- Knee hyperextension

Solutions

- Shift socket anteriorly in relation to the foot
- Attach stiffening bumpers. See **Stiffening Bumpers** section below for installation instructions

Heel too hard

Symptoms

- Rapid knee flexion, instability
- Heel to toe progression too rapid
- Lack of energy return sensation

Solutions

- Shift socket posteriorly in relation to the foot
- Verify appropriate foot module category

Foot module too stiff

Symptoms

- Flat spot in rollover motion at slow cadences

Solutions

- Consider a lower category foot module

Foot module too soft

Symptoms

- Clicking noise at *initial contact*
- Excessive toe deflection with high impact activity

Solutions

- Consider a higher category foot module

Stiffening Bumpers

Stiffening bumpers are included to adjust the heel stiffness during *loading response*. The bumpers may be temporarily attached between the heel lever and the keel using the pre-applied adhesive in the location indicated on the bumper package to increase heel stiffness one category. If the heel stiffness is too stiff, move the bumper posteriorly; still too soft, move it anteriorly. For permanent placement, clean off the pre-applied adhesive with Acetone, and adhere bumpers using Super Glue (cyanoacrylate).

Spectra™ Socks

One Spectra™ sock is provided to minimize noise and to protect the foot shell/glass composite components. The Spectra™ sock should be placed over the keel and the sole plate before donning the foot shell. Spectra™ socks must be replaced at intervals appropriate to the user's activity level. Failure to inspect and replace the Spectra™ socks may prematurely wear the foot module, and will void the warranty.

Foot Shell

When removing or installing the foot shell, use the Foot Shell Removal Tool (ACC-00-10300-00) to prevent damage to the foot module.

Maverick

Maximum user weight: 166 kg (365 lbs.)
Available sizes: 23cm-30cm

**Instructions for use
Heel Bumper**

Test bumper placement

Temporarily attach at recommended location using pre-applied adhesive. Suggested placement. If necessary, adjust anterior-posterior.

Permanently bond:

Remove temporary adhesive on bumpers using a cloth with acetone. Do not soak bumpers in acetone. Bond bumpers in pre-selected position using cyanoacrylate (Super Glue).

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Advena Limited, Tower Business Centre, 2nd Flr., Tower